




**FRESENIUS
KABI**

caring for life

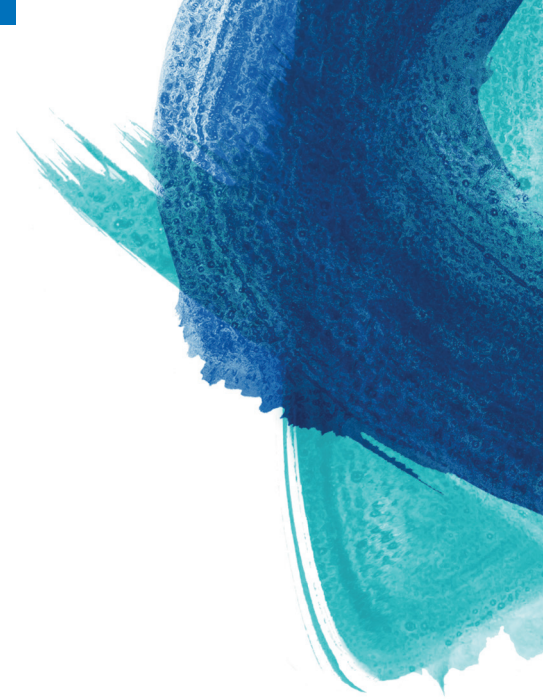

Donor Tube Stripper

Used to strip undiluted blood from donor tubing into blood collection container.

Ergonomic design

Light weight machined aluminum

Durable electroless nickel coating with satin finish

Adjustable roller gap

Donor Tube Stripper

The Fresenius Kabi Donor Tube Stripper is crafted with durable, lightweight machined aluminum. The rollers are supported with reinforced pins for reliable performance.

The Donor Tube Stripper includes a 2mm hex key and 0.9mm gap gauge for adjusting the roller gap.



Description	Product Code	Weight	Closed Dimensions
Donor Tube Stripper	6R4452	0.5 lb (227 grams)	6.75" L x 1.75" W x 1" H (17.1 cm x 4.5 cm x 2.5 cm)

For more information please contact us at ProductSupport@fresenius-kabi.com or 800-937-5060.

Refer to the Instructions for Use for a complete list of warnings and precautions associated with the use of this product.

Trademarks referred to are property of their respective owners.
Copyright © 2022 Fresenius Kabi AG. All rights reserved.